

Cordless Water Flosser Instruction Manual

Model: FJ3001A

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

WARNINGS

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger. The water flosser is only to be charged with the input DC 5V charging adaptor.
- Water flosser can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure they do not play with the product.
- Product shall be used only as intended, following the instructions outlined in this manual.
- The water flosser is a personal care appliance and is not intended for use on multiple patients in a dental practice.
- Never operate this product if it has a damaged cord, plug or charger. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the product, including power cord, away from heat surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- Do not attempt to remove battery. Attempts to remove battery will void the warranty.
- Do not clean the water flosser in the dishwasher.
- Keep the charger out of the reach of children at all times.

DANGER-WATER

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not place or store the charger where it can fall or be pulled into a bathtub, sink or toilet.
- Do not place or drop the charger into water or other liquid.
- Do not reach for a product that has fallen into water or other liquid. Unplug immediately.
- Unplug charger before cleaning and ensure that it is dry before plugging it into an electrical outlet.

DANGER-VOLTAGE

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug the charger into an outlet with a voltage other than specified on the bottom of the charger.
- This product is designed to operate within a range of DC 5 volts.
- Voltage converters do not guarantee compatibility.
- If the plug does not fit easily into the outlet, discontinue use.
- Never wrap the cord tightly around the appliance. Never allow the power supply cord to be pulled, twisted, or severely bent. Untwist the cord before use.

MEDICAL WARNINGS



- If you have had oral surgery or gum surgery within the past two months, consult your dentist before using this product.
- Consult your dentist if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- Consult your physician prior to using the water flosser if you have any medical concerns.
- To avoid injury, point nozzle and spray only between teeth. Do not directly spray gums or throat. It is recommended that you turn the unit off when switching from upper teeth to lower.

HOUSEHOLD USE ONLY

- This product has no user serviceable parts.
- Store the charging unit in a dry and safe place when it is not in use.
- Untwist the power cord before use.
- Do not dismantle or exchange/manipulate any part of the water flosser if it malfunctions or fails. This will void the warranty and you may run the risk of an electric shock.
- If the supply cord is damaged, it must be replaced by the manufacturer in order to avoid a hazard.
- If the product is dropped, the tip should be replaced before the next use even if no damage is visible.

SAVE THESE INSTRUCTIONS

PRODUCT DESCRIPTION

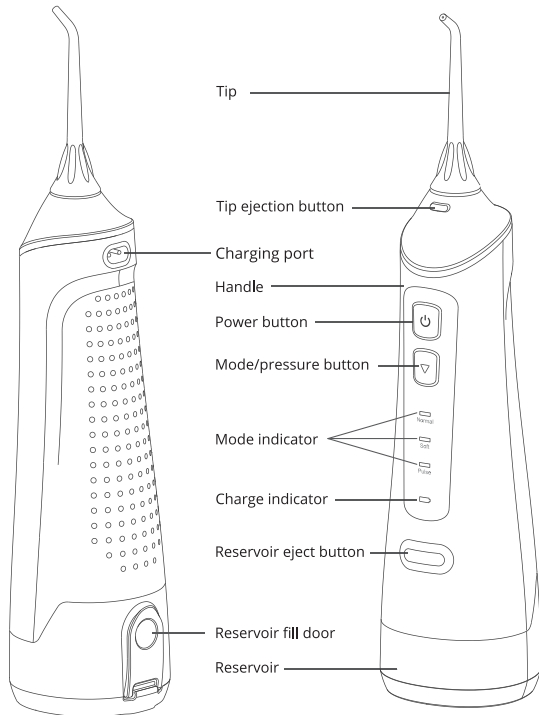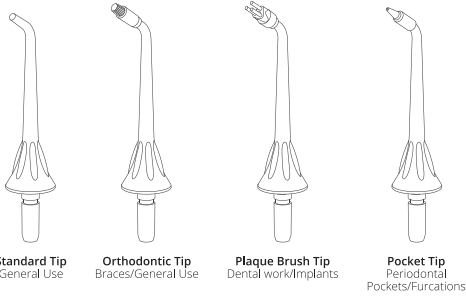

CLEANING MODES

The water flosser automatically starts in the default mode.

NORMAL

- Standard mode for superior teeth cleaning.

SOFT

- Gently cleans your teeth or gums and other sensitive areas.

PULSE

- Pulsates for gentle gum stimulation.

HOW TO USE

Charging the water flosser

Please fully charge prior to first use.

• Plug the power cord into the USB socket and plug the cord on the back charging port of the unit (Fig.1). The battery charge indicator flashes to indicate that the unit is charging. When the unit is fully charged, the battery charge indicator stops flashing and remains on.

Note:

- When the batteries are in need of charging, the charge indicator will flash off after 5 seconds.
- A complete recharge will take 4-5 hours. A full charge will last for approximately 40 times of typical use.
- During charging, the unit handle will get warm, but this is normal.
- The operating duration will become shorter as battery gets older. If the unit operates only for a short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact customer service for support.

Inserting and Removing Tips

- Insert tip into the center of the knob at the top of the water flosser handle (Fig.2).
- Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle (Fig.3).

CAUTION: Do not press tip eject button while unit is running.

Filling the Reservoir

- Remove the charger cable from the unit prior to filling reservoir and product use. Lift flip top on reservoir and fill the reservoir with lukewarm water (Fig.4).

Adjusting the Pressure Setting







You can adjust the pressure setting on your water flosser by pressing the mode button on the handle. Start with the lowest setting and increase to your liking (Fig.5).

Recommended Technique

- Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON (Fig.6).
- For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.
- Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

When Finished

- Turn the unit off. Empty any liquid left in the reservoir.

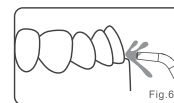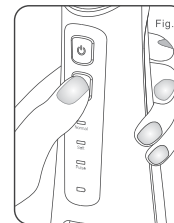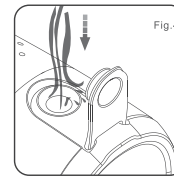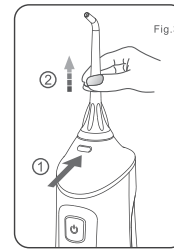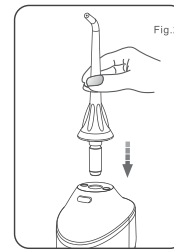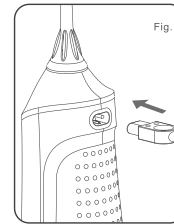

TIP USAGE

Standard Tip - General Use

The standard tip is for everyday use. It cleans deep between teeth and below the gumline to help prevent gum disease and improve your oral health.


To use the standard tip, start with the back teeth and follow the gumline, pausing briefly between the teeth and aiming the tip at the gumline at a 90-degree angle. Hold the tip slightly away from the gums to allow the water flow to remove plaque and debris.



Orthodontic Tip - Braces/General Use

The orthodontic tip has a tapered brush on the end to help remove plaque from braces and orthodontics. It also helps flush out bacteria and food debris from around teeth and under the gums.

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

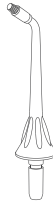

Plaque Brush Tip - Dental work/implants

The plaque brush tip removes plaque around all types of dental work and hard-to-reach areas. It has 3 thin tufts of bristles that gently access stubborn plaque around dental implants, crowns, bridges, and veneers.

Place the tip close to the teeth so that the bristles are gently touching the teeth. As you water floss, gently glide the tip along the gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.



Pocket Tip - Periodontal Pockets/Furcations

The pocket tip is specifically designed to deliver water or antibacterial solutions deep into periodontal pockets.

To use the Pocket tip, **set the unit to the lowest pressure setting.** Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.




CLEANING

- Do not clean handle or replaceable tip in the dishwasher or a microwave.
- Never immerse unit in water. Before cleaning, unplug from the electrical outlet.
- The water reservoir is removable for easy cleaning and is top-rack dishwasher safe. It is recommended that the reservoir is removed and cleaned weekly.
- Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.


Cleaning the tip

- Rinse the tip after each use.
- Remove the tip from the handle and rinse the tip connection at least once a week with warm water (Fig.7).

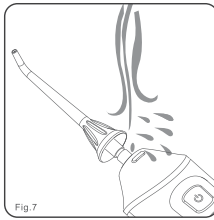

Storage

- If you are not going to use the water flosser for an extended period of time, unplug the charger, clean it and store it in a cool and dry place away from direct sunlight.
- If desired, the reservoir may be removed from the power handle by pressing the reservoir release latch to the unlock position, and then sliding the reservoir down towards the base of the unit (Fig.8).

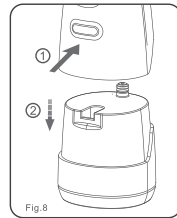

TROUBLESHOOTING

If you experience an issue, we are here to help!

| PROBLEM | POSSIBLE CAUSE | SOLUTION |
|--|---|--|
| Water does not come out. | There is no water in the reservoir or unit is tilted too much. | Ensure the reservoir is full and the straw in the reservoir is sitting in the water. |
| The water flosser doesn't work. | You have just purchased the product or it has not been used for more than 6 months. | Charge the handle for at least 5 hours. |
| The water flosser operates only for a short period of time even after being charged. | The battery has reached the end of its life cycle. | Contact customer service for support. |
| | Charging time was insufficient. | Charge the handle for at least 5 hours. |

LIMITED TWO-YEAR WARRANTY



Smile Brilliant warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. When you purchase from the Smile Brilliant website, your purchase receipt is automatically saved. You will need your order number to reference any warranty claims. For a warranty claim, it may be requested that you return the defective product at your expense. Smile Brilliant will replace at its discretion any part of the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Smile Brilliant. This limited warranty excludes accessories or consumable parts such as tips, etc.

Use with improper voltage will damage product and void warranty.

FCC STATEMENT


This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference
- (2) This device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment.

UNIT AND BATTERY DISPOSAL

At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. Always remove the battery before discarding.

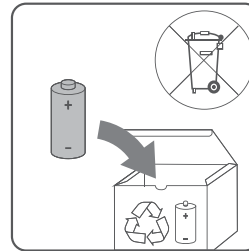

CONTACT



**Monday - Friday
9 a.m - 5p.m CST**



Call us at **855.944.8361**



Email us at **Support@SmileBrilliant.com**



Live chat with us at **www.SmileBrilliant.com**